



# Fully Automated

T.D.M, D.O.A & Immunoassay



# Compact Dedicated Analyser at affordable Cost



*Quality and Precision!*

# FEATURES

90 Tests / Hour  
20 Reagent Positions  
Disposable Cuvettes  
Autoloading and Disposal  
User Interfaced Built in LCD  
Barcoded Reagents  
Cooled Reagents

# SPECIFICATIONS

Type:	Desktop fully automated random access analyser	
Sample type:	Serum, plasma and urine	
Array type:	Endpoint and rate arrays, 1 and 2 point	
Calibration type:	Factor, linear, 2 point and multipoint	
Sample input system:	Removable tray with sample tube holder on a turntable Barcode sample ID reading with optional handy barcode reader	
Capacity:	Up to 10 samples (handling of STAT available)	
Pipettor:	Number of pipettors One for sample and reagent	
Type:	Micro-pipette with liquid level detector	
Rinsing:	Inside and outside with purified water	
Sampling pump:	Type:	Micro syringe
	Volume:	2 to 35 $\mu$ l (increments of 0.1 $\mu$ l)
Reagent system:	Removable tray with reagent bottle holder on a turntable Barcode reagent ID reading with optional handy barcode reader	
Capacity:	20 reagent positions on a tray (ten 50ml and ten 20ml bottle positions)	
Inventory:	Calculation of remaining reagent volume available	
Cooling:	Cooling with Peltier element (8 to 15° C)	
Reagent pump:	Type:	Micro syringe
	Volume:	20 to 350 $\mu$ l (increments of 1 $\mu$ l)
Cuvette system:	Number of cuvettes:	24 cuvettes on reaction line
	Type:	Disposable with automatic loading system
	Material:	Resin
	Volume for measurement:	Minimum 1-20 $\mu$ l Maximum 350 $\mu$ l in total
Reaction system:	Type:	Direct heating system
	Temperature:	37° C
Dector:	Method:	Direct measurement of absorbance in cuvette Bi-chromalic and Mono-chromatic
	Wavelength:	340 to 800, 8 discrete wavelengths by filter
Light source:	Halogen tungsten lamp	
Stirring system:	Stick type rotating stirrer	
Dimension:	620 mm (W) x 555 mm (D) x 570 mm (H)	
Weight:	65kg	
Option:	Hand held barcode reader	

## • HEADQUARTERS IRELAND



Audit Diagnostics.  
Business & Technology Park,  
Carrigtwohill,  
Co.Cork,  
Ireland.  
Tel +353 21 4533652  
Fax +353 21 4533653  
Email: [info@auditdiagnostics.ie](mailto:info@auditdiagnostics.ie)  
Website: [www.auditdiagnostics.com](http://www.auditdiagnostics.com)

## • CHINA



Audit Diagnostics China  
Building 16A-1406  
Fangchengyuan 1 Zone,  
Fengtai District  
Beijing 100078, China.  
Tel / Fax: +86-1391 056 5056  
Email: [minfenglin@auditdiagnostics.ie](mailto:minfenglin@auditdiagnostics.ie)

## • INDIA



Audit Diagnostics India Pvt Ltd.  
Chanderwati Vashisht Kunj  
1st Floor  
247, Jain Nagar  
Karala, Delhi-110081  
India.  
Tel:+91(11)25474733/65164533  
Fax:+91(11)25968849  
E.Mail: [vashisht@auditdiagnosticsindia.com](mailto:vashisht@auditdiagnosticsindia.com)  
web: [www.auditdiagnosticsindia.com](http://www.auditdiagnosticsindia.com)

## • TURKEY




Audit Diagnostics Turkey  
Ruzgarli Cad. Eser Han  
Kat. 3 No. 14/175  
Ulus, Ankara, Turkey.  
Tel: +90 312 309 48 61  
Fax: +90 312 309 48 62  
Email: [auditkemal@gmail.com](mailto:auditkemal@gmail.com)

## • UK



Audit Diagnostics UK Ltd.  
P.O.Box 679,  
Wivelsfield Green,  
Haywards Heath,  
West Sussex, RH 16 9BY,  
UK.  
Tel / Fax + 44 1444 471 494  
Email: [selwyn\\_morgan@msn.com](mailto:selwyn_morgan@msn.com)

Distributor

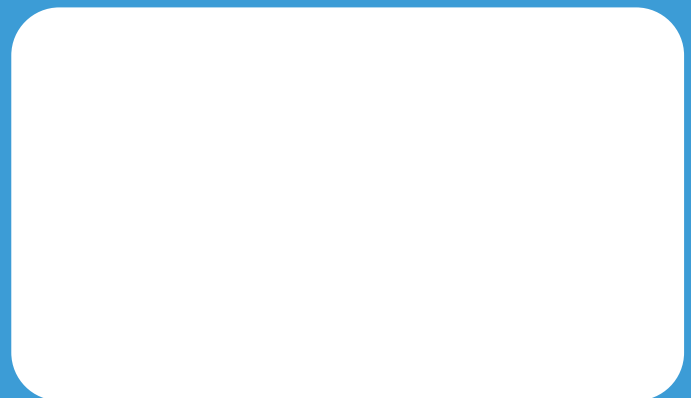

To Order or For Further Information  
Contact



**Audit Diagnostics**